




NRG LIGHT

NrG-AURORA Stand & Pulley Guide

Setup and safety guide for the included rolling stand and door-mount pulley system.

NrG-AURORA

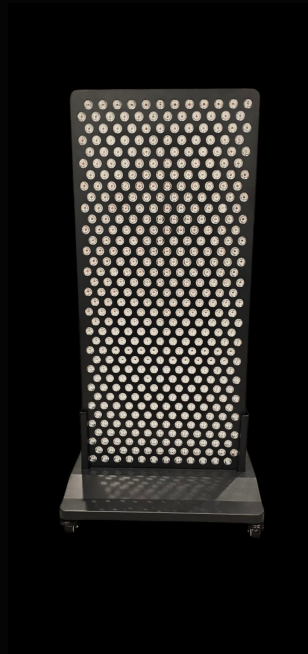

Revision 1.0 • August 2026

ROLLING STAND

Before mounting the panel

- Assemble the stand only on a level floor using the hardware supplied with the stand.
- Lock all casters before lifting or attaching the Aurora.
- Use two people when attaching or removing the panel so one person can support the panel while the other secures the mounting hardware.
- Tighten all supplied fasteners evenly and confirm the panel cannot shift in the bracket before use.
- Route the power cable so it remains clear of wheels, adjustment points, and walking paths.

USING THE ROLLING STAND

- Keep casters locked during every session and whenever adjusting panel position.
- Turn the panel off before moving the stand. Unlock the casters only when you are ready to roll the unit.
- Push from the stand structure at a controlled walking speed; do not pull the panel by its enclosure or cable.
- Avoid thresholds, steep transitions, or uneven flooring that could tip the stand.

Use only on a suitable door and with the supplied hardware

- Inspect the door, hinges, frame, hooks, pulley, cord, and attachment hardware before each setup.
- Use only a structurally sound door that closes securely. Do not use a damaged, loose, or unsuitable door.
- Seat the supplied over-door hardware fully and make sure the door cannot swing unexpectedly during use.
- Attach the panel only to the intended mounting points and keep the panel supported until the pulley system is carrying the load securely.
- Raise and lower the panel slowly. Keep hands, hair, clothing, and cords clear of the pulley and pinch points.
- Before standing or lying beneath the panel, test the full setup at the intended height and confirm that every connection is secure.

Important: Mounting hardware can vary by production run. The supplied hardware diagram and the physical hardware included with your unit take precedence. If anything does not match this guide, stop and contact NrG Light before assembly.

SUPPORT

Questions about stand or pulley setup?

Contact info@nrglight.com. For panel controls, use the shared NrG Light Panel User Manual.